<b>D</b>	<b>Pureo Derma Peel</b> <b>Mikrodermabrasion</b> Gebrauchsanweisung .....	2	<b>I</b>	<b>Pureo Derma Peel</b> <b>Microdermoabrasione</b> Istruzioni per l'uso .....	34
<b>GB</b>	<b>Pureo Derma Peel</b> <b>microdermabrasion</b> Instructions for use .....	10	<b>TR</b>	<b>Pureo Derma Peel</b> <b>Mikrodermabrazyon</b> Kullanım Kılavuzu .....	42
<b>F</b>	<b>Pureo Derma Peel</b> <b>Microdermabrasion</b> Mode d'emploi .....	18	<b>RUS</b>	<b>Pureo Derma Peel прибор по</b> <b>уходу за кожей лица</b> Инструкция по применению .....	50
<b>E</b>	<b>Pureo Derma Peel</b> <b>Microdermoabrasión</b> Instrucciones de uso .....	26	<b>PL</b>	<b>Urządzenie do mikrodermabrazji</b> <b>Pureo Derma Peel</b> Instrukcja obsługi .....	58

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89231 Neu-Ulm

Germany

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

## Table of contents

1. Getting to know your instrument .....	11
2. Signs and symbols .....	11
3. Intended use .....	11
4. Warnings and safety notes .....	12
5. Unit description .....	13
6. Use .....	14
6.1 Carrying out a skin test .....	14
6.2 Using the device .....	14
6.3 Application areas .....	15
7. Maintenance and cleaning.....	16
8. Frequently asked questions.....	16
9. Disposal.....	16
10. Technical data.....	16
11. Replacement parts and wearing parts .....	17

## Included in delivery

- 1x FC 100 base station
- 1x hand-held device
- 1x precision attachment
- 1x coarse attachment
- 1x large attachment
- 20x single-use filter
- 1x mains adapter
- These instructions for use

## WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device is not suitable for persons under the age of 18.
- This device may be used by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

## Dear customer,

Thank you for choosing one of our products. Our brand stands for well-engineered, quality made products; driven by innovation and safety. Our product range encompasses a wealth of products for health and well-being in the heat, gentle therapy, blood pressure/diagnosis, weight, massage, beauty and air sectors.

With kind regards,  
Your Beurer team.

## 1. Getting to know your instrument

The Pureo Derma Peel FC100 enables professional facial peels with microdermabrasion technology. A high-quality sapphire coating is used to gently remove excess skin cells and the state-of-the-art vacuum technology noticeably improves skin texture and enhances circulation. The treatment also simulates cell renewal for a finer and smoother complexion.

## 2. Signs and symbols

The following symbols appear in these instructions for use and on the device:

	<b>WARNING</b> Warning instruction indicating a risk of injury or damage to health
	<b>IMPORTANT</b> Safety note indicating possible damage to the device/accessory
	<b>DANGER</b> This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!
	<b>Note</b> Note on important information
	Observe the instructions for use.
	The device is double-protected.
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Manufacturer
	Only use indoors.

## 3. Intended use

The device is only to be used on the face (apart from the eye area and mouth) of adults. Do not use the device on animals!

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

## 4. Warnings and safety notes

### **WARNING**

- Danger of suffocation by packaging material. Keep packaging material away from children.
- Before use, inspect the device and accessories for visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.
- Keep the mains adapter away from water.

### **DANGER**

This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!

### **WARNING**

Do **NOT** use the device

- on cracked skin
- on open wounds
- in the eye area or on the mouth
- on any part of the body other than the face (except for eye area and mouth)
- if you have a skin disease (e.g. herpes, vascular diseases or skin growths)
- if you have skin cancer
- if you are suffering from skin irritation
- immediately after sunbathing or if you have sunburn
- if you are taking medication that contains steroids
- if you are taking blood thinners
- if you have multiple sclerosis, diabetes or an autoimmune disorder

### **WARNING**

For this reason, the device must only be operated as follows only with the supplied power supply unit at the mains voltage indicated on the power supply unit.


### **IMPORTANT**

Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. In the event of a fault or damage, have the device repaired by a certified specialist workshop.

### **WARNING**

- Only use the device in dry rooms (e.g. never in the bath, sauna).
- Never reach for a device that has been submerged in water. Disconnect the mains adapter immediately.
- Do not drop the device or subject it to any impacts.
- Do not clamp or pinch the mains cable.
- Ensure that the mains cable and mains adapter do not come into contact with water, steam or other liquids.

## 5. Unit description

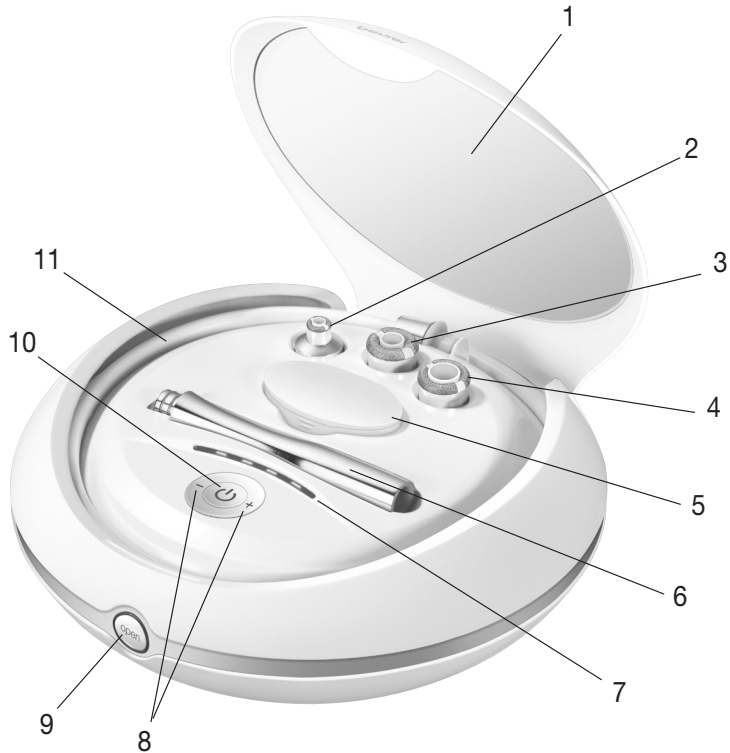1. Mirror	7. Intensity level display (level 1-5)
2. Precision attachment with fine sapphire coating, suitable for precise treatment of small areas (e.g. in nose area)	8. +/- buttons for setting the intensity level
3. Large attachment with fine sapphire coating, suitable for sensitive skin	9. Open button
4. Coarse attachment with coarse sapphire coating, suitable for intensive treatment of problem areas	10. ON/OFF button
5. Storage box with single-use filters	11. Tube
6. Hand-held device	

## 6. Use

### Note

Before you start using the device, make sure that your face is clean, dry and free of any traces of products such as powder, cream or make-up.

### 6.1 Carrying out a skin test

Carry out a skin test 24 hours before using the device. Test the device at the lowest intensity level on small, inconspicuous areas of your neck. If after 24 hours the tested area is red or painful, do not use the device.

### WARNING

- Do not use the device more than once a week.
- Immediately after use, the treated skin may be slightly reddened and in very rare cases there may be small spots of bleeding under the skin.
- For one day before and one day after use, avoid activities that cause you to perspire, such as saunas, tanning beds and sun exposure.


### Note

The device can also be used twice a week for non-sensitive problem skin (with a break of at least 48 hours between use).


### 6.2 Using the device

1. Connect the mains adapter to the device and plug the adapter into the socket.


2. Press the **Open button** to open the device.

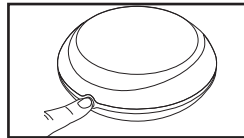
3. Insert a single-use filter into the end of the tube.

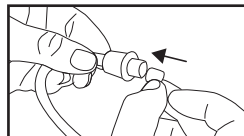4. Connect the end of the tube to the hand-held device.

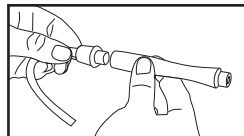
5. Place the attachment you wish to use on the tip of the hand-held device.

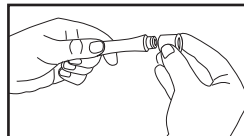6. Now press the **ON/OFF button** to switch on the device. The intensity level display starts to light up. You are now at the lowest intensity level.

7. Start the treatment at the lowest intensity level initially. Place the hand-held device on your face and move it over the skin. Make sure that you do not move the hand-held device over one part of the face more than once. Further details on treatment can be found in Section "6.3 Application areas".

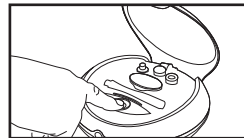
If the lowest intensity level is not strong enough for you, you can change this by pressing the **-/+ buttons**.

### **i** Note

The pressure you apply when using the device should be comfortable and not painful. Some slight reddening should be expected after use.

8. When you have finished using the device, press the **ON/OFF button** to switch it off and unplug the mains adapter from the socket.

9. Change the single-use filter after every use. To do this, disconnect the hand-held device from the tube. Use tweezers to remove the single-use filter.
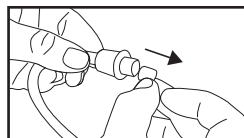
## **6.3 Application areas**

### **i** Note

**Ensure that the part of your face you are treating is taut when you move the hand-held device over the skin.** You can do this by pulling the skin taut to the side or downwards with your free hand.

#### **Forehead:**

Move the hand-held device slowly and evenly outwards from the centre of the forehead. Make sure that you do not treat any part of your skin more than once.

#### **Temples:**

Move the hand-held device outwards towards your hairline. Repeat this to the left and right of the nose. Ensure that the device is at least 1 cm away from your eyes.

#### **Nose:**

Move the hand-held device downwards from the centre at the top of the bridge of your nose. Repeat this to the left and right of the nose.

#### **Cheeks:**

Move the hand-held device outwards across your cheeks.

#### **Around the mouth:**

Move the hand-held device along your upper lip from the centre to the left and right. Move the hand-held device upwards from your bottom lip to the left and right. Make sure that you do not touch your lips.

## 7. Maintenance and cleaning

- Never submerge the device in water or other liquids.
- Disinfect the attachments with a standard disinfectant after each use.
- Change the single-use filter after each use.
- Use a dry cloth to clean the device of any dirt. Do not use any abrasive cleaning products.
- Protect the device from knocks, damp, dust, chemicals, marked temperature fluctuations, and nearby sources of heat (ovens, heaters).

## 8. Frequently asked questions

### Is the Pureo Derma Peel FC100 suitable for all skin types?

The Pureo Derma Peel FC100 is suitable for all skin types. It can be used on oily, dry and combination skin. It is also suitable for all skin colours.

### Is the technology in the Pureo Derma Peel FC100 safe?

The Pureo Derma Peel FC100 is based on safe microdermabrasion technology. This only removes excess skin cells in the very top layer of the skin and does not affect the other skin layers.

### Is using the Pureo Derma Peel FC100 painful?

No, use of the Pureo Derma Peel FC100 should be gentle and pain-free. On very sensitive skin, users may experience slight burning in rare cases. If, however, you do experience pain when using the device, reduce the intensity.

### How long does a treatment with the Pureo Derma Peel FC100 take?

The length of a treatment can vary depending on the area being treated. A full facial treatment takes approx. 20 minutes.

### How often should I use the Pureo Derma Peel FC100?

We recommend using the Pureo Derma Peel FC100 once a week and waiting at least 48 hours between treatments. You can also use the device twice a week on non-sensitive problem skin. A test must be carried out before use. Protect your skin from the sun after use.

### What precautions need to be taken after treatment?

Before exposing yourself to the sun, use sun protection for the first three days after a treatment. For the first three days after a treatment DO NOT carry out any mechanical/chemical exfoliation treatment, electrolysis or facial hair removal and DO NOT use hair removal products. You can wear make-up immediately after treatment.

## 9. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point.

Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local authorities responsible for waste disposal.

## 10. Technical data

Mains part voltage supply	Input: 100-240V AC, 50/60 Hz Output: 12V, 600mA
Dimensions	Ø 21,5 x 8,8 cm
Weight	~869 g

## 11. Replacement parts and wearing parts

Replacement parts and wearing parts are available from the corresponding listed service address under the stated material number.

<b>Designation</b>	<b>Item number and/or order number</b>
20 single-use filters	163.546
Precision attachment	163.557
Coarse attachment	163.556
Large attachment	163.555

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Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany